

Features

- UL® 1069 Listed
- Provides for the following functions from a Hill-Rom® SideCom® Unit:
 - ◇ Call-placed indication
 - ◇ Monitor indication
 - ◇ Hands-free two-way communication
 - ◇ Nurse Call
- Supervised
- Ramp lock plug-in connectors


Description

The RY351B SideCom® adapter is designed for use with IR311E, IR312E, IR319B and IR320B patient stations, and a Hill-Rom® SideCom® unit. When used with the above patient stations and the Hill-Rom® SideCom® unit, the RY351B adapter will provide all the necessary circuitry for a Hill-Rom® SideCom® unit to control all functions normally controlled by a pillow speaker. This includes, but is not limited to, nurse call and patient communications. Note: There is one RY351B SideCom® adapter required per SideCom® unit.

Listings



UL® 1069 listed; UL® File Number E71599.

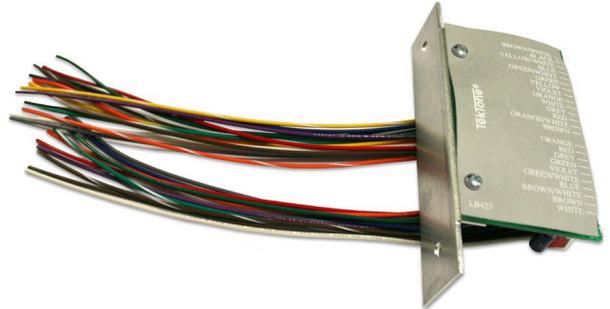

RY351B Hill-Rom® SideCom® Adapter

Specifications

Dimensions:	Length:	4.25"	(108 mm)
	Width:	2.00"	(50.8 mm)
	Height:	1.25"	(31.8 mm)
Construction:	Circuit board mounted on a rugged aluminum mounting bracket.		
Connections:	Reliable ramp lock plug-ins.		
Housing:	Single gang ring, or single gang ring and double gang box, or exact equivalent.		

Required Components

Tek-CARE®300II Nurse Call System
JK401 Dummy plug for IR319B or IR320B patient stations