



Senior Living

Wireless Pendant

Designed by residents and their caregivers, the senior living pendant lets residents live their most active life without compromising safety or style.

TekTone[®]
COMMUNICATIONS FOR TODAY ... AND TOMORROW.

Tek-CARE[®] 500

 **INOVONICS**[®]



Designed by Residents and Caregivers

The senior living pendant provides residents with the comfort of knowing help is just a button press away, while offering the advanced functionality and reliability caregivers and communities depend on.

Resident Focused

60%

Smaller and lighter than our previous pendant



Residents have high expectations for their pendants. The senior living pendant delivers the fashion and functionality they demand.

Small, Light and Comfortable

60% smaller and lighter than our previous pendant, the new pendant is discreet and easy to conceal under clothing. The contoured shape and polished finish make it comfortable to wear.

IP67 Waterproof Rating

Sealed housing and water tight battery compartment allow residents to safely shower or bathe while wearing their pendant.

Advanced Functionality

Bright lights signal alarm and low-battery conditions. Vibration feature confirms alarm activation.

Enhanced Wearability

Both neck lanyard and wristband options are specifically tailored for comfort. An optional belt loop kit is also available.





Caregiver Approved

Caregivers need a pendant that's dependable and easy to service. The senior living pendant allows them to focus on residents instead of their devices.

Proven Reliability

Functions over the proven EchoStream® wireless network and is certified for installation in UL 2560 listed emergency call systems.

Robust Battery Life


Unique power management specifically designed for senior living usage patterns provides one year of battery life at 1800 activations (five/day).

Innovative Alarm Clearance

Choice of magnetic alarm clearance card or button press pattern eliminates reset button. Bright light signals successful alarm clearance.

Waterproof Twist-Lock Battery Door

Battery changes are quick and easy with no screws. The available alarm clearance card also unlocks the battery door.

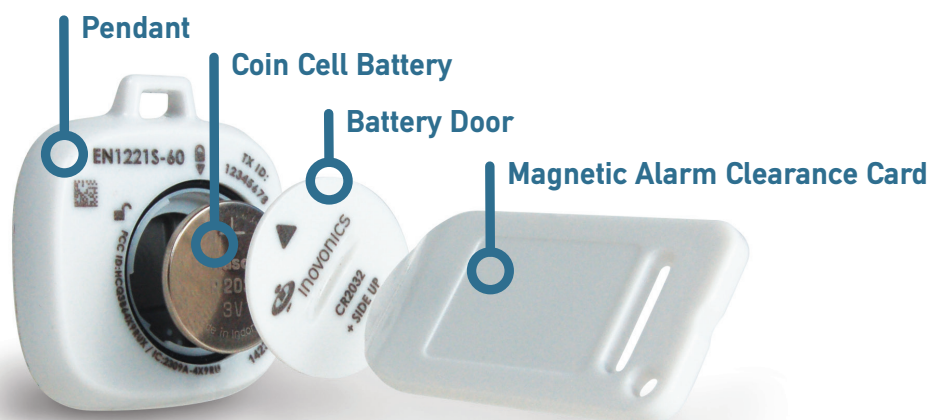

Packed full of functionality

An amazing amount of new features in one incredibly small package

Large Alarm Light

Bright and continuous to confirm activation

Vibration Confirmation

Immediate assurance of successful activation

High Contrast Button

Easy to see, with raised dome to locate by touch

Waterproof Housing

Sealed housing achieves IP67 waterproof rating for use in shower or bath

Twist-Lock Access

The replaceable coin cell battery is accessed with a simple twist

Attachment Options

Versatile keyring and retractable cord slot attachment points

Low Battery Light

Provides five day warning to residents and staff

Smart Power Management

State of the art battery system provides continuous power and one year battery life

Simple Alarm Clearance

Alarms can be cleared by placing the magnetic card over the pendant button

Advanced Multi-Function Lights



Alarm Activation

When activated, pendant vibrates briefly and light flashes red until alarm is cleared.



Alarm Clearance

Light flashes blue when magnetic card is used or button pattern is pressed.



Low Battery

Light blinks yellow to notify residents and staff.

Pendant Accessories



Magnetic alarm clearance card



Replacement neck lanyard



Replacement wristband



Belt loop kit for neck lanyard pendant



Replacement battery door