

### Features

- 24" touchscreen monitor for call annunciation and interaction
- Can be placed anywhere in the facility with LAN connectivity
- Touchscreen monitor for easy interaction
- Displays calls, icons, maps, station information, resident information, and much more
- Calls can be filtered and zoned
- Can operate as a non-interactive display if only a call display is desired
- Optional supervision
- Desk or wall mountable. Wall mount uses a standard VESA mounting system. Wall mounts provided by contractor


### Description

The touchscreen NC403TS Tek-CARE Monitor is a purpose-built monitor for the Tek-CARE nurse call system. It is designed to be a cost-effective and attractive option for displaying calls throughout the facility where the full functionality of a master station is not required.

The NC403TS Tek-CARE Monitor is supplied as a complete unit which includes all the necessary hardware for connection to the nurse call system. The NC403TS is connected to the Tek-CARE system via the facility's local area network (LAN).

The NC403TS Tek-CARE Monitor offers many functions used by facility staff on a daily basis, such as call display and the ability to view details about a call, resident information, status display, and more. In addition to a call list view, a map view and icon view are also available if desired. The NC403TS Tek-CARE Monitor can also modify its own annunciation settings and display system faults if desired. Password protection of advanced features is also available.

The NC403TS Tek-CARE Monitor can be zoned to display calls from desired zones to improve staff awareness and efficiency. Up to 254 NC403TS Tek-CARE Monitors can be installed on a single Tek-CARE system.

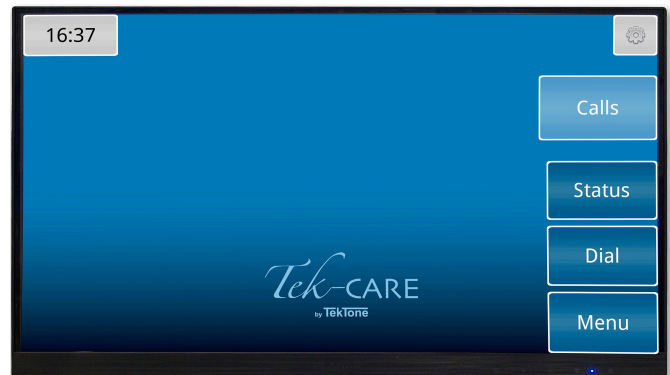

NC403TS Tek-CARE Monitor

### Specifications

- Dimensions: 14"H × 21.3"W × 1.9"D  
(356 mm × 541 mm × 48 mm)
- Connection: Ethernet connectivity to facility LAN
- Display: 24" touchscreen display

### Required Components

NC475 Tek-CARE Appliance Server